


LikoScale™ 200

Accessory

User Guide*



Prod. No. 3156225

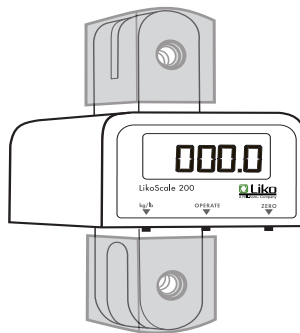

LikoScale 200 accessory is a convenient, easy-to-use scale that is intended for weighing patients in Liko mobile or stationary lifts. Together with Liko's lifts and lifting accessories, LikoScale 200 accessory enables safe and easy weighing.

The weight is shown on the scale's digital display. By selecting unit mode the scale can display either kg or lbs.

The scale should only be used in combination with approved lifts and accessories in accordance with the instructions in this guide.

LikoScale 200 accessory can be used with the following Liko lifts:

Viking™ mobile lift (S, M, L, XL)

- require Adapter 12 mm, Prod. No. 2016504.

Uno™ mobile lift (all models)

- Uno models with serial numbers from 31825 require Adapter 12 mm, Prod. No. 2016504.
- Uno models with serial numbers up to 31824 require Adapter 8 mm, Prod. No. 2016502.

Golvo™ mobile lift (all models)

Likorall™ overhead lift (all models)

Multirall™ overhead lift (all models)

Liko™ MasterLift overhead lift (all models)

As the lifting height is limited when LikoScale accessory is used in combination with LikoLight™ mobile lift or Viking™ XS mobile lift, we do not recommend these alternatives.

* This user guide is a complement to the complete documentation issued by the manufacturer of the scale. The intention with this document is to make it easier for our customers to use LikoScale in combination with Liko's products. For detailed information regarding the scale we refer you to the manufacturer's Operation Guide (attached).

Before assembly, make sure you have the following:

- Tool kit, 2 screws, 2 locking nuts
- 2 sleeves
- Battery
- User Guide from Liko
- Operating and Service Manual from manufacturer
- Declaration of Calibration

Installing/replacing the battery


LikoScale 200 accessory is delivered with a 9 V battery with a battery life of about 3000 readings.

Installation of battery:

1. Open the door on the back of the scale.
2. Install the battery.
3. Close the door.

When the battery needs to be replaced, a small battery symbol will appear in the upper left corner of the display.

NOTE! Old batteries should be left at the nearest recycling station.

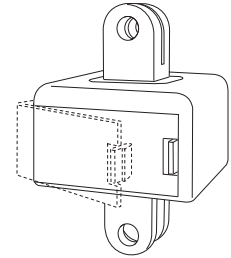

Assembly Instructions

LikoScale 200 accessory is delivered with two sleeves. The attachment type you have above and below the scale, and the type of sling bar you have determines if both sleeves must be used. Follow the assembly instructions for the appropriate lift, attachment point, and sling bar.

Viking™ mobile lift, Uno™ mobile lift: Fixed assembly with adapter

1. Remove the lifting accessory from the lift.
2. Mount the appropriate adapter (see previous page) on to the lift in accordance with the assembly instruction.
3. Attachment point under the scale:

Alternative 1 with Universal sling bar



Place one sleeve on the scale's lower attachment point.

Mount the lifting accessory under the scale, using bolt and locking nut provided.

Alternative 2 with Standard sling bar

Mount the lifting accessory under the scale, using bolt and locking nut provided.

4. Place one sleeve on the scale's upper attachment point. Mount the scale in its upper attachment to the adapter, using bolt and locking nut provided.

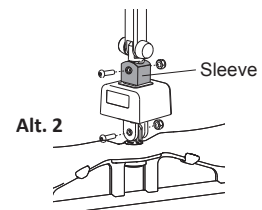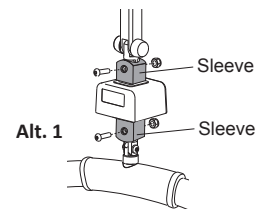


△ When disassembling the scale, the adapter must be removed from the lift before a lifting accessory is reassembled on to the lift.

Viking™ mobile lift, Uno™ mobile lift: Quick-release assembly

Liko lifts equipped with (Q-Link) quick fastener require the Adapter Kit LikoScale, Prod. No. 3156232.

1. Remove the lifting accessory from the lift.
2. Place one sleeve on the scale's upper attachment point.
3. Mount the Adapter Kit LikoScale in accordance with the assembly instruction.
4. Mount the scale on to the lift.
5. Mount the lifting accessory under the scale.

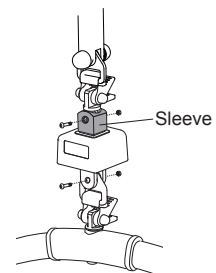

Likorall™ overhead lift, MasterLift overhead lift, Golvo™ mobile lift: Fixed assembly

1. Remove the lifting accessory from the lift.
2. Attachment point under the scale:

Alternative 1 with Universal sling bar



Place one sleeve on the scale's lower attachment point.

Mount the lifting accessory under the scale, using bolt and locking nut provided.

Alternative 2 with Standard sling bar

Mount the lifting accessory under the scale, using bolt and locking nut provided.


3. Mount the scale on to the lift, using bolt and locking nut provided.

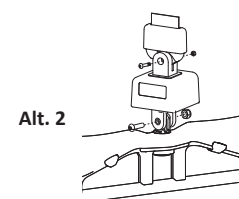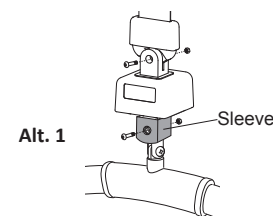

Multirall™ overhead lift, Likorall™ overhead lift, MasterLift overhead lift, Golvo™ mobile lift: Quick-release assembly

Liko lifts equipped with (Q-Link) quick fastener require the Adapter Kit LikoScale, Prod. No. 3156232.

1. Remove the lifting accessory from the lift
2. Place one sleeve on the scale's upper attachment point.
3. Mount the Adapter Kit LikoScale in accordance with the assembly instruction.
4. Mount the scale on to the lift.
5. Mount the lifting accessory under the scale.

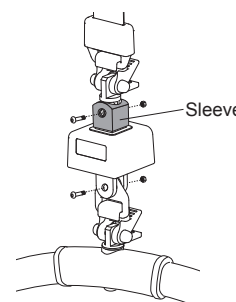

Operation

The normal climate range for operation of LikoScale 200 accessory is between +10° C (50° F) to +40° C (104° F) at a relative humidity not exceeding 95%.

The scale must not be stored at temperatures below -20°C (-4°F) or above +50°C (+122°F).

When the scale is stored for more than 6 months, the battery should be removed. Follow instructions for Installing/replacing the Battery.

Before weighing, ensure the scale has been installed in accordance with the Assembly Instructions on pages 2-3.

⚠ Before lifting, check that the lifting accessory is hanging vertically and can move freely.

Setting of unit mode (kg/lbs.)

Press the (kg/lbs) - button on the scale to display either kg or lbs.

Weighing with LikoScale 200 accessory

The scale can be reset to zero with or without the sling. If the scale is reset to zero without the sling attached to the sling bar, the weight of the patient will include the weight of the sling. Therefore the weight of the sling must be subtracted to obtain the correct weight of the patient.

Resetting the scale with the sling is described below:

1. Connect the sling to the sling bar.
2. Start the scale by pressing the "OPERATE" button.
3. Reset the scale by pressing the "ZERO" button.
When the scale is ready, the value "0.0" will appear on the display.
4. Remove the sling from the sling bar and apply the sling to the patient in accordance with the instruction guide for the sling and lift. Connect the sling to the sling bar.
5. Raise the sling bar so the patient rises slightly. Never lift a patient higher than necessary.
6. When the lifting motion has stabilized, press the "OPERATE" button. The patient's weight can now be read on the display.
The scale turns off automatically after about 60 seconds. For further readings, press the "OPERATE" button.
7. Lower the patient carefully and remove the sling from the sling bar.

Calibration

For calibration, see manufacturer's Operating and Service Manual or contact Liko.

Technical Data

Display resolution:	0.1 lb / 0.1 kg
Accuracy:	0.1 % +/- Digit of reading
Max. load:	200 kg (440 lbs.)
Min. load:	10 kg (22 lbs.)
Reset:	One button reset ("ZERO")
Power supply:	9 V Battery
Display:	Digital, liquid crystal (LCD)
Size:	76 x 76 x 102 mm (3 x 3 x 4 in.) (W x L x H)
Weight:	0.6 kg (1.3 lbs.)
Temperature range:	+10° C to +40° C (+50° F to +104° F)
Humidity range:	0% - 95%



Distributor:

Liko AB
Nedre vägen 100
SE-975 92 Luleå, Sweden

Manufacturer:



SR Instruments, Inc
600 Young Street
Tonawanda, NY 14150, USA

Enhancing outcomes for
patients and their caregivers:

Hill-Rom