


# Hill-Rom® Trauma Stretcher

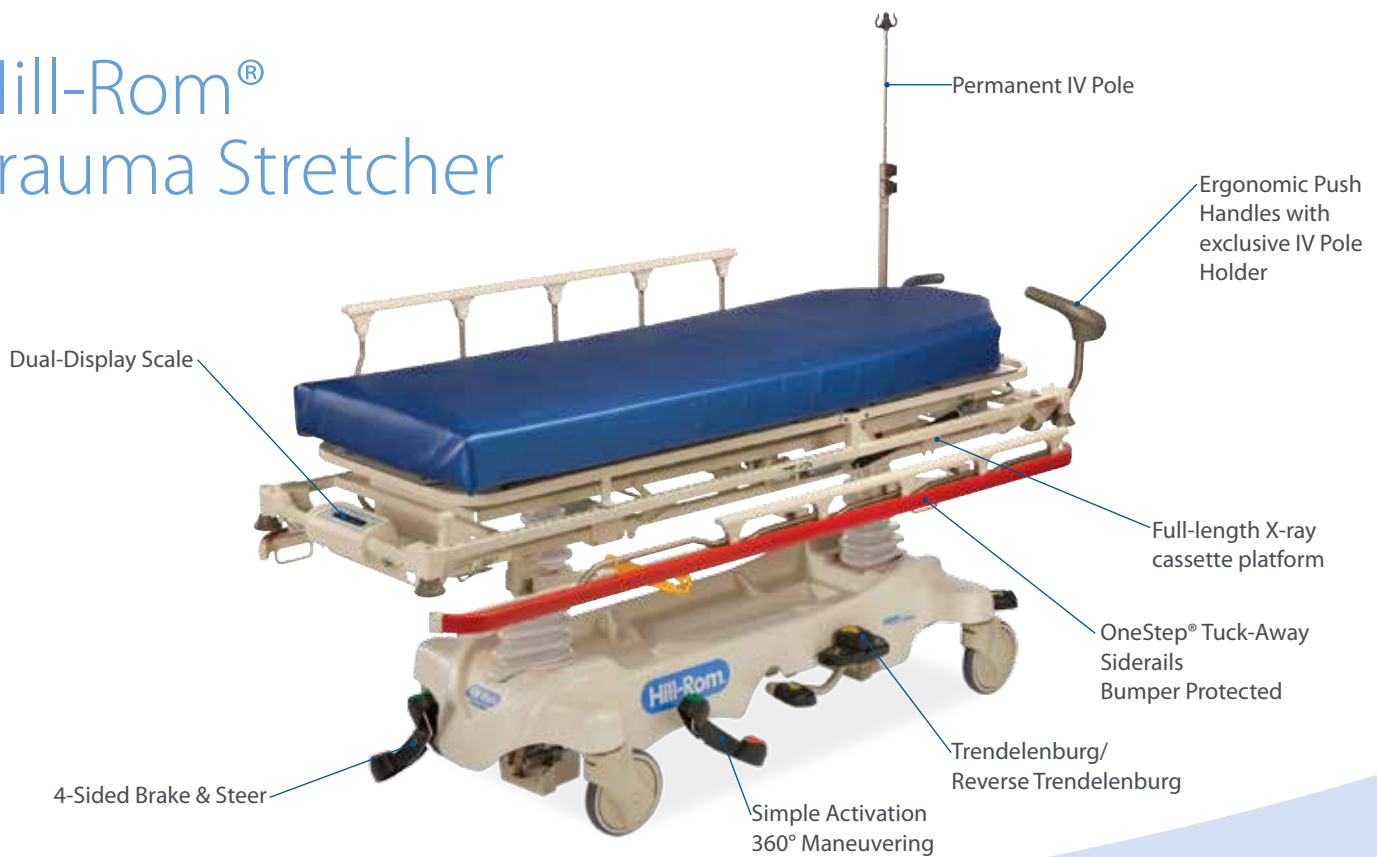

## Go Beyond Mobility™

Our Trauma Stretcher provides versatility to help you care for your high acuity patient. With the built-in x-ray cassette platform, capture quality images right on this stretcher, without moving your patient.

A rush to the operating room is not just about running down hallways and turning corners. It's about the minute details that contribute to a safe and comfortable experience for both you and your patient.

Our stretchers **provide more than mobility** for you and your patient:

- A safe and comfortable patient experience
- Simple and safe maneuvering
- Durability to withstand your demands

## Key benefits at a glance:

### A Safe and Comfortable Patient Experience

- Help prevent pressure injuries with your choice of pressure redistribution surfaces.
- Protect your patient's hands with bumper protected siderails.
- Provide a gentle patient ride with spring-loaded Steering Plus™ system that absorbs shock of driving over bumps and thresholds.

### Simple and Safe Maneuvering

- Accommodate standard-size X-ray cassettes with full-length cassette deck.
- 360° maneuvering with Steering Plus™ system
- Zero Transfer Gap with OneStep® Tuck-Away Siderails
- Obtain an immediate patient weight with scale that memorizes last weight.

### Durability to Withstand Your Demands

- 700 lb weight capacity
- Compatibility with most lifts and overbed tables
- Easy to clean smooth deck and siderail design
- Power washable

# Go Beyond Mobility™

A safe and comfortable experience for both you and your patient starts on a Hill-Rom® stretcher.



Cassette Locking Lever



AccuMax® VPC Surface



Steering Plus™ System



Dual-Display Scale



Ergonomic Push Handles with Exclusive IV Pole Holder



Active Hand Brake

A SAFE & COMFORTABLE PATIENT EXPERIENCE	STANDARD FEATURES	OPTIONAL FEATURES
26" width	X	
Bumper Protected Siderails	X	
Trendelenburg/Reverse Trendelenburg	X	
Six Transport Strap Locations	X	
PRESSURE REDISTRIBUTION MATTRESS OPTIONS		
3" Comfortline® foam mattress	X	
4" or 5" Comfortline® foam mattress		X
5" AccuMax® VPC Surface (Elastimax or Hypolex)		X
SIMPLE & SAFE MANEUVERING		
Steering Plus™ System	X	
OneStep® Tuck-Away Siderails	X	
Brake/Steer – head and foot	X	
Brake/Steer – all four sides		X
Dual-Sided hi/lo Hydraulic Control Pedals	X	
3-Sided hi/lo Hydraulic Control Pedals		X
Ergonomic Push Handles		X
Ergonomic Push Handles with Integrated IV Pole Transport Holder		X
Active Hand Brake		X
3/4" Standard IV Pole Sockets	X	
Removable IV Pole – 2 Stage		X
Permanent IV Pole – 3 Stage (40 lb weight capacity)		X
Integrated Dual-Display Scale		X
Radiolucent Deck	X	
DURABILITY TO WITHSTAND YOUR DEMANDS		
700 lb weight capacity	X	
Full Perimeter Bumping	X	
Smooth Deck Design	X	
Decals Options	X	
Bumper Color Options – choice of neutral, blue, red, teal or purple	X	

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2017 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.  
197942 rev 2 30-JAN-2017 ENG – US

For further information about this product or a service, please contact your local Hill-Rom representative or visit our webpage:

USA 800-445-3730  
Canada 800-267-2337

[www.hill-rom.com](http://www.hill-rom.com)

## Technical Specifications

Maximum weight limit	700 lbs (317.5 kg)
Low position (floor to top of sleep deck)	24.25" (61.6 cm)
High position (floor to top of sleep deck)	37.5" (95.3 cm)
Overall length	83" (210.8 cm)
Overall width (siderails up)	32" (81.3 cm)
Overall width (siderails stored)	29.25" (74.3 cm)
Siderail length	48" (121.9 cm)
Siderail height above sleep deck	10.7" (21.2 cm)
Comfortline® mattress size	26" (66 cm) x 75" (190.5 cm)
AccuMax® VPC mattress size	26" (66 cm) x 76" (193 cm)
Maximum head elevation	90°
Maximum Trend/Reverse Trend	18°
Floor to base clearance	3.5" (8.9 cm)
Caster size	8" (20.3 cm)

## Warranty

Frames & welds for life of product

Three-year parts

Two-year mattress

One-year labor

Accessory Options — Please confirm accessory availability for model of stretcher ordered

Accessory IV pole	Liquid oxygen tank holder
Ankle stirrups	Oxygen tank holder
Armboard	Paper roll dispenser
Chart holder	Patient tray
Convertible footboard	Restraint straps
Detachable footboard	Siderail covers
Foot extender	Stretcher pillow
Infusion support system	Transport straps
IV transporter	Utility tray
Lateral X-ray cassette holder	

Enhancing outcomes for patients and their caregivers:

**Hill-Rom**