




Welch Allyn  
SureTemp® Plus 692 Electronic Thermometer

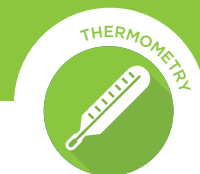

# The SureTemp Plus 692 Thermometer— with the speed and accuracy you require, and the safety and security you demand.

1

## Easy:

Intuitive user interface, easy-to-read LCD screen and flashing icons, helps select and identify temperature measurements and modes.

2

## Safe:

Designed with interchangeable oral/axillary (blue) and rectal (red) probe wells to address patient safety and help reduce the risk of cross-contamination.

3

## Fast and Accurate:

With proven accuracy, the SureTemp 692 thermometer captures quick readings—approximately 4-6 second oral, 10-13 second pediatric axillary and rectal readings and 12-15 second adult axillary readings.

4

## Secure:


To top it all off, our security features make sure the unit stays where it is supposed to be. A location ID on the display and other innovative features help prevent theft.

- Easy-to-use, intuitive operation
- Interchangeable oral/axillary and rectal probe wells help reduce the risk of cross-contamination



# SureTemp Plus 692—setting the standard of care

ID field for electronic labeling of thermometer

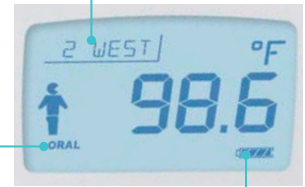

User-selectable icon and/or words for identification of thermometer mode

Battery life indicator

Interchangeable probe with push-button probe cover ejection

Removable probe well for easy cleaning

C/F button for toggling between Fahrenheit and Celsius

Prominent recall button for display of last temperature taken

Ergonomic design with rubberized grips

Convenient storage for 25 probe covers

Large, easy-to-read, backlit LCD

Mode button for selecting oral, pediatric axillary or adult axillary

Pulse timer button/monitor mode

**NEW** feature makes it easier to switch to Monitor Mode for measuring a patient's thermal steady state (approx. 3 to 5 minutes)

Interchangeable rectal probe and well reduce risk of cross-contamination



The **SureTemp Plus 692** with **standard wall mount** has storage for an additional 50 probe covers.



# The SureTemp Plus 692—setting the standard of care even higher

## SureTemp Plus 692 Product Specifications

<b>Calibration Accuracy</b>	±0.2 °F (±0.1 °C) (Monitor Mode)	
<b>Temperature Predict Time</b>	<b>Oral:</b> <b>Adult Axillary (18 years and older):</b> <b>Pediatric Axillary (17 years and younger):</b> <b>Rectal:</b>	Approximately 4-6 secs Approximately 12-15 secs Approximately 10-13 secs Approximately 10-13 secs
<b>Unit Dimensions</b>	8.46" x 3.18" x 2.43" (215 mm x 81 mm x 62 mm)	
<b>Weight</b>	12.6 oz (357 g)	
<b>Ambient Temperature Range</b>	50.0 °F (10 °C) to 104 °F (40.0 °C)	
<b>Patient Temperature Range</b>	80.0 °F (26.7 °C) to 110 °F (43.3 °C)	
<b>Power Source</b>	(3) 1.5Vdc AA batteries	
<b>Battery Life</b>	Approximately 5,000 temperatures	

## Ordering Information:

### THERMOMETERS

- 01692-200** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Security System with ID Location Field, Oral Probe and Oral Probe Well with 4ft cord
- 01692-201** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Security System with ID Location Field, Rectal Probe and Rectal Probe Well with 4ft cord
- 01692-300** SureTemp Plus 692 Electronic Thermometer, One per room, Wall Mount, Oral Probe and Oral Probe Well with 9ft cord (must be locked to the wall)
- 01692-301** SureTemp Plus 692 Electronic Thermometer, One Per Room, Wall Mount, Rectal Probe and Rectal Probe Well with 9ft cord (must be locked to the wall)
- 01692-700** SureTemp Plus 692 Electronic Thermometer, Wall Mount, Oral Probe and Oral Probe Well with 9ft cord, with Rolling Stand and Holder for Spare Probe & Well

### ACCESSORIES

- 21326-0000** Wall Holder
- 21330-0000** Holder for Spare Probe and Well
- 06138-000** Calibration Key Assembly, 690/692
- 02892-000** Probe and Well Kit, 4ft, Rectal
- 02892-100** Probe and Well Kit, 9ft, Rectal
- 02893-000** Probe and Well Kit, 4ft, Oral
- 02893-100** Probe and Well Kit, 9ft, Oral
- 02891-0000** Probe Well, Blue
- 02891-1000** Probe Well, Red
- 08273-000** Thermometer Stand with Holder for Spare Probe and Well
- 21333-0000** Bracket, M600 Stand Upgrade, with Holder for Spare Probe and Well
- 98025-0000** Extended Warranty M692

### REPLACEMENT PROBE COVERS

- 05031-101** Case of 1,000 Probe Covers

Contact your Welch Allyn representative or visit [www.welchallyn.com](http://www.welchallyn.com) to learn more.



Welch Allyn, Inc.  
 4341 State Street Road  
 Skaneateles Falls, NY 13153 USA  
 (p) 800.535.6663 (f) 315.685.3361