

Hill-Rom® Electric Stretcher



Go Beyond Mobility™

Our Electric Stretcher provides the comfort of a hospital bed for those unexpected long stays. With positioning controls integrated into the siderails, patients are in control of their own comfort.

A rush to the operating room is not just about running down hallways and turning corners. It's about the minute details that contribute to a safe and comfortable experience for both you and your patient.


Our stretchers **provide more than mobility** for you and your patient:

- A safe and comfortable patient experience
- Simple and safe maneuvering
- Durability to withstand your demands

A Safe and Comfortable Patient Experience

Patients are spending more time on stretchers, and as a result, the patient experience, including pressure injury prevention, has become increasingly important.

- **Help prevent pressure injuries** through your choice of pressure redistribution surfaces.
- **Customize comfort** with electric patient positioning controls.
- **Prevent patient migration** with electric knee articulation.
- **Protect your patient's hands** with bumper protected siderails.
- **Provide a gentle patient ride** with spring-loaded Steering Plus™ system that absorbs shock of driving over bumps and thresholds.
- **Allow for safe patient ingress and egress** with low height.
- **Provide a comfortable option** for patients to recover with Recovery Position.

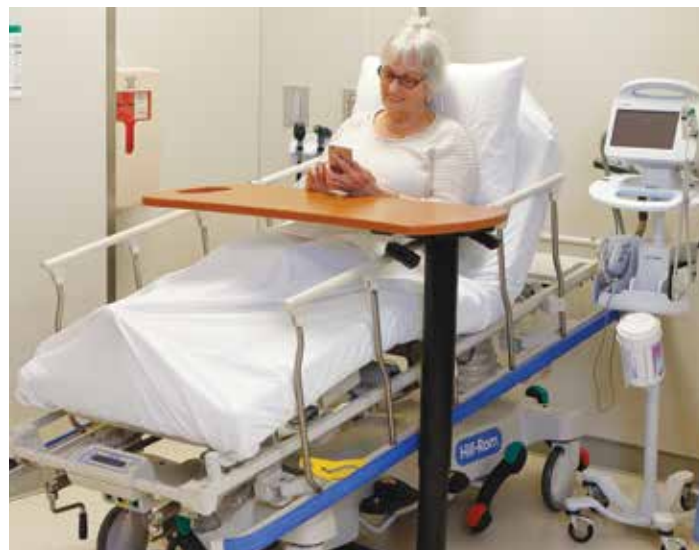

Recovery position





Simple and Safe Maneuvering

Whether you're making crisp turns down narrow halls, getting the stretcher rolling from a stop, or nestling it into a crowded elevator, our stretchers are ideal in any setting. Hill-Rom Stretchers are easy to maneuver to help reduce your risk of injuries.

- **Simple activation** requires less force to activate than Stryker's Big Wheel.
- **360° maneuvering** with Steering Plus™ system grips the floor and reduces your turning radius.
- **Ease the burden off your back** with electric positioning controls.
- **Zero Transfer Gap** — OneStep® Tuck-Away Siderails for seamless patient transfers.
- **Obtain an immediate weight without moving your patient** with scale that memorizes last weight.

Durability to Withstand your Demands

At Hill-Rom, we understand that because caregivers have to do more with less, having a reliable and durable product is more important than ever.

- **700 lb weight capacity** for all stretchers
- **Compatibility** with lifts and overbed tables
- **Easy to clean** smooth deck and siderail design

Go Beyond Mobility™

A safe and comfortable experience for both you and your patient starts on a Hill-Rom® stretcher.



Patient Controls



Caregiver Control Panel



AccuMax® VPC Surface



Dual-Display Scale



Steering Plus™ System



Active Hand Brake

| A SAFE & COMFORTABLE PATIENT EXPERIENCE | STANDARD FEATURES | OPTIONAL FEATURES |
|---|-------------------|-------------------|
| 30" width | X | |
| Electric Fowler Articulation | X | |
| Electric Knee Articulation | X | |
| Bumper Protected Siderails | X | |
| Trendelenburg/Reverse Trendelenburg | X | |
| Six Transport Strap Locations | X | |
| PRESSURE REDISTRIBUTION MATTRESS OPTIONS | | |
| 4" Comfortline® foam mattress | X | |
| 5" Comfortline® foam mattress | | X |
| 5" AccuMax® VPC Surface (Elastimax or Hypolex) | | X |
| SIMPLE & SAFE MANEUVERING | | |
| Steering Plus™ System | X | |
| OneStep® Tuck-Away Siderails | X | |
| Caregiver Electric Fowler Control | X | |
| Caregiver Electric Knee Control | X | |
| Caregiver Fowler & Knee Lockout Controls | X | |
| Brake/Steer – head and foot | X | |
| Brake/Steer – all four sides | | X |
| Dual-Sided hi/lo Hydraulic Control Pedals | X | |
| Ergonomic Push Handles | | X |
| Ergonomic Push Handles with Integrated IV Pole Transport Holder | | X |
| Active Hand Brake | | X |
| 3/4" Standard IV Pole Sockets | X | |
| Removable IV Pole – 2 Stage | | X |
| Permanent IV Pole – 3 Stage (40 lb weight capacity) | | X |
| Integrated Dual-Display Scale | | X |
| DURABILITY TO WITHSTAND YOUR DEMANDS | | |
| 700 lb weight capacity | X | |
| Full Perimeter Bumpering | X | |
| Smooth Deck Design | X | |
| Decals Options | X | |
| Bumper Color Options – choice of neutral, blue, red, teal or purple | X | |

Hill-Rom reserves the right to make changes without notice in design, specifications and models. The only warranty Hill-Rom makes is the express written warranty extended on the sale or rental of its products.

©2017 Hill-Rom Services, Inc. ALL RIGHTS RESERVED.
197892 rev 2 10-FEB-2017 ENG – US

Technical Specifications

| | |
|---|---|
| Maximum weight limit | 700 lbs (317.5 kg) |
| Low position (floor to top of sleep deck) | 22.5" (57.2 cm) |
| High position (floor to top of sleep deck) | 34.25" (87 cm) |
| Overall length | 83" (210.8 cm) |
| Overall width (siderails up) | 36" (91.44 cm) |
| Overall width (siderails stored) | 30.5" (77.5 cm) |
| Siderail length | 48" (121.9 cm) |
| Siderail height above sleep deck | 14.5" (36.8 cm) |
| Comfortline® mattress size | 26" (66 cm) x 75" (190.5 cm) or 30" (76.2 cm) x 75" (190.5 cm) |
| AccuMax® VPC mattress size | 30" (76.2 cm) x 75" (190.5 cm) |
| Maximum head elevation | 90° |
| Maximum Trend/Reverse Trend | 18° |
| Floor to base clearance | 3.5" (8.9 cm) |
| Caster size | 8" (20.3 cm) |

Warranty

Frames & welds for life of product

Three-year parts

Two-year mattress

One-year labor

Accessory Options — Please confirm accessory availability for model of stretcher ordered

| | |
|-------------------------------|---------------------------|
| Accessory IV pole | Liquid oxygen tank holder |
| Ankle stirrups | Oxygen tank holder |
| Armboard | Paper roll dispenser |
| Chart holder | Patient tray |
| Convertible footboard | Restraint straps |
| Detachable footboard | Siderail covers |
| Foot extender | Stretcher pillow |
| Infusion support system | Transport straps |
| IV transporter | Utility tray |
| Lateral X-ray cassette holder | |

For further information about this product or a service, please contact your local Hill-Rom representative or visit our webpage:

USA 800-445-3730
Canada 800-267-2337

www.hill-rom.com

Enhancing outcomes for
patients and their caregivers:

Hill-Rom